relaxmyway

CELLUMAX

< relaxmyway

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
Copyright © 2016 CASADA International GmbH.

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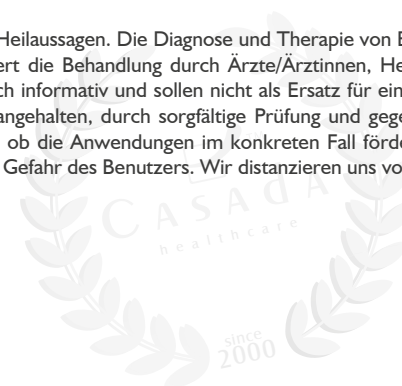< relaxmyway

Since the company was founded in the year 2000, the name CASADA has stood for high grade products in the area of wellness and fitness.

Today, CASADA – products are available worldwide. What distinguishes CASADA products is their continuous development, unmistakable design and their absolute focus on improving quality of life as well as their reliable and excellent quality.

Congratulations!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Cellumax.

Table of Contents

10	Safety instructions
10	Start up
11	Technical data
11	Equipment / Functions
12	Application
12	Troubleshooting
13	Storage
13	Conformity declaration
13	Warranty conditions

Safety Instructions

Before using this massage device, please read the instructions carefully to ensure perfect function and optimal effect. Please retain this instruction manual in a safe place!

- The massage device complies with the recognised principles of technology and the current security regulations.
- The equipment is maintenance free. Any incidental repairs may only be executed by an authorised expert. Inappropriate usage as well as unauthorised repairs are prohibited due to reasons of safety and will lead to the loss of warranty.
- Avoid touching the plug with wet hands.
- Please avoid the equipment coming in contact with water, high temperatures as well as direct exposure to the sun.
- In case of a longer period of non-usage please unplug the equipment from the power supply to avoid risk of short circuit.
- Do not use damaged plugs, cables or loose sockets. If the plug or the cables are damaged, they have to be replaced by the manufacturer, the service representative or qualified personnel. In case of malfunction please immediately interrupt the connection with the power supply. Misappropriation or wrongful usage will void all liability for damages.
- Children should only use the massage device under appropriate supervision.
- In order to prevent excessive stress on muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Never use pointy or sharp objects.
- Any massage – also massage by foot – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: recent injuries, thrombotic illnesses, inflammations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of afflictions and ailments.
- If you depend on electrical aids such as a pace maker, please consult your physician prior to usage in any case.

Start Up

CONNECTING THE CABLE

Please check prior the power plug and connect it with the power network. Now turn on the device by the control wheel.

Technical Data Cellumax

Measure: 19 x 10,5 x 17,5 cm
Weight: 1 kg
Power source: AC 220 – 240 V, 50 Hz
Rated power: 25 Watt

Certificates:

Equipment / Functions

CelluMax – massage with a deep effect. Two different, easily interchangeable massage heads active the nervous and circulatory system in two different ways. The blood circulation of the massaged area is activated.

Easy handling: thanks to the control wheel the speed of the CelluMax can be adjusted without any level stages and make a more individual massage possible. The integrated heating function intensifies the deep effect. CelluMax can be used together with e.g. cellulite creme to attain better results. The cover is removable and washable. The special curved handle with its anti-slip coating was designed to ensure easy and safe utilization.

[1] CONTROL WHEEL ON/OFF/INTENSITY

[2] FLAT HEAD

This top part massages gently and is suitable for sensitive areas such as the face (please use with cover since little hairs to get entangled)

[3] NUBBLED HEAD

The shape of this top part pleases with an invigorating message: the nervous and circulatory system are being activated, the blood circulation supported. Especially suitable to use with cellulite cream (please remove cover before).

Application

CelluMax massages the connective tissue and supports blood circulation and the elimination of toxins. It affects the body fat, warms it up and drains it. By using the CelluMax you bring your body in shape, massage tired and tense body parts and overcome cellulite.

Nevertheless, you have to be conscious about your weight. Overweight will cause an increased storage of fat in your fat cells. Therefore a balanced diet is very important. Exercises for abs, butt and thighs will tighten your muscles and help to burn fat. With only 20 minutes a day you can already achieve remarkable results.

By massaging the certain affected areas blood circulation and metabolism are activated. After showering, legs and buttocks should be showered with cold water. Consider also that through nicotine blood vessels narrow and that the metabolism is hindered.

Troubleshooting

ERROR

The massage device can not operate even though it has been connected to the power supply.

ADJUSTMENT

- Check where the plug is well connected and working well.
- The massage device could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn OFF the power and allow the device to cool off for 40 to 50 minutes.

Storage

1. Recommended storage: dry, dust-free. Please roll up the wire.
2. Do not store next to a heater or open flame. Do not expose to direct sun.

Conformity Declaration

Corresponds to the European guidelines EMC 2014/30/EU, LVD 2014/35/EU.

Warranty Conditions

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt. The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorized dealer.

The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medical treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
info@casada.com
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright © 2017 Casada International GmbH. All rights reserved.