TM
CASADA
healthcare

< relaxmyway

TAPPYMED II

< relaxmyway

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< relaxmyway

Since being established in 2000, the name CASADA stands for high quality products on the health-spa and fitness market. Today, CASADA-equipment is available in 37 countries world-wide.

Constant further development, an unmistakable design and a focused orientation towards a better quality of life are what characterize CASADA-products, as much as their unbeatable and reliable quality.

Congratulations!

With the purchase of this massage device you have shown awareness of your health.

In order for you to enjoy the advantages of your massage device in the long term, we ask you to read and adhere to the instruction manual carefully.

We wish you a lot of fun with your personal Tappymed II.

Table of Contents

Safety instructions	12
Equipment / Functions.....	13
Technical Data	13
Massage technique	14
Applications	14
Operating instructions.....	15
Troubleshooting.....	16
Care instructions	16
Storage	16
Conformity declaration	16
Warranty condition.....	17

Safety Instructions

Please carefully read the following instructions before using this massage device to ensure its proper function and optimal effect. Please retain these instructions for use!

- This massage device complies with the recognised principles of technology and current safety regulations.
- This equipment is maintenance free. Any necessary repairs must be performed by an authorised professional. Improper use and unauthorised repairs are prohibited for safety reasons and will void the warranty.
- Never touch the power plug with wet hands.
- Please avoid the equipment coming into contact with water, high temperatures and direct sunlight.
- Be sure to unplug the equipment from the power point during extended periods of non-use to eliminate the risk of a short circuit.
- Do not use damaged plugs, cords or loose sockets. Damaged plug or cords must be replaced by the manufacturer, service representative or qualified personnel. In the event of malfunction unplug the equipment immediately. We assume no liability for damages resulting from wrongful or improper use.
- In order to prevent excessive stress on muscles and nerves, the recommended massage time of 30 minutes should not be exceeded.
- Never use pointy or sharp objects.
- Any massage – also massage by hand – has to be avoided during pregnancy and/or in case of one or more complaints in the area of the massage: recent injuries, thrombotic illnesses, inflammations and swelling of any kind as well as cancer. Consultation of your physician is recommended prior to any treatment of afflictions and ailments.
- Do not use the device while driving.
- If you depend on electrical aids such as a pace maker, please consult your physician prior to usage in any case.

Equipment / Functions

- Relieves tension in the back, waist and entire body.
- Relaxes taut muscles and tense nerves, especially after workouts.
- Alleviates aches and pains in the arms, neck and shoulders.
- Stylish, anti-slide grip design makes the massager easy to handle.
- Specially designed streamline design enables easy kneading to any part of the body.
- Adjustable speed for effective, powerful massage to stimulate nerves and muscles.
- Focused massage helps to relieve muscle aches, induce relaxation and revitalize the body.
- Pounding effect relaxes tight muscles after doing exercise while easing tendons and relieving congestion.
- Heating function
- Four shapes of massage caps for different massage styles.
- Red and Blue LED light help to distinguish the different functions of the unit easily: Red LED indicates that the Vibrating and Heat function is working. Blue LED indicates that the Vibrating function is working.

Technical Data Tappymed II

Measurements device:	10 x 40 x 9 cm
Weight:	990 g
Power source:	200~240V /50-60 Hz
Rated Power:	20 Watt
Speed:	high 2600 revolutions per minute low 2200 revolutions per minute
Operating time automatic:	15 minutes
Certificates:	   

Massage technique

TAPPING

During the Tapping short knocking movements are performed with the edge of the hand, the flat hand or the fist. They enhance the circulation of the skin and change the tension in the muscles. If Tapping is performed at the level of the lungs, it will relieve mucus build-up in the lungs.

Applications

- Relieves pain and muscle tightness in the back.
- Relaxes tired and stressed muscles, especially after sporty activity.
- Relieves tension and cramps in Shoulders and neck area.
- Prevention through targeted leg massage or reflexology.

5 DIFFERENT WAYS TO USE THE TAPPYMED II:

- [1] Without a head: This is like tapping with a ball, and can be used to stimulate especially the lower part of the fist, the palm of the hand, or the arms.
- [2] The “brush” like head: This gives a softer massage and is well suited for the joints of the body like: knees, elbows, shoulders, ankles, feet. As the massage is softer, you may find this head to be a good and comfortable alternative to be used for any parts of your body.
- [3] The “nipple” like head: It focuses the tapping motion to just one point. You can use it particularly on those parts of the body which are likely to have tense muscles. E.g. in the thigh region. Provides a deeper massage. It does also work perfectly well as a reflexology massage for the sole of your feet, as you can pinpoint the different reflexology points using the “nipple”

[1] Without a head [2] brush

[3] nipple

[4] 3-Point


[5] bar

Applications

- [4] With the “3-Point” head: Works in a similar way as the “nipple” head, but focuses on 3 points for a more rigorous massage. Very effective for e.g. the calf, using the massage under pressure. It can also be applied softly to the joints of your fingers.
- [5] With the “bar” like flat head: Very effective massage for long muscular regions like the back and shoulder area, as well as the hamstrings.

Operating instructions

- [1] Insert power plug into the outlet.
- [2] Select massage speed.
- [3] Set the speed switch to OFF if you wish to stop kneading.
- [4] Press once on HEAT button to activate the heat function, press again to stop the heat function. While the red light is working, the heat function is ON.
- [5] Do not massage longer than 15 minutes at a time. Using the massager longer than 30 minutes could over-stimulate your muscles and nerves, leading to negative effects.
- [6] Please remove the power plug from the outlet.

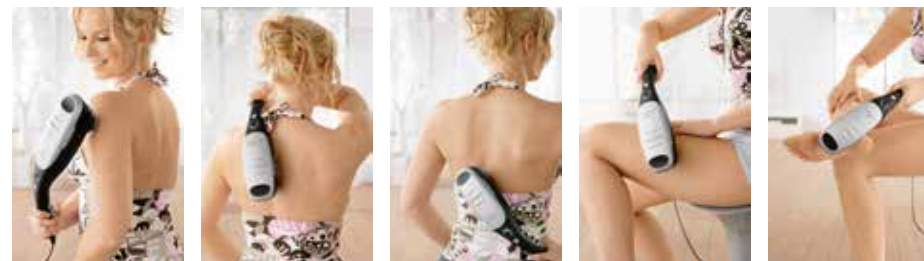Troubleshooting

ERROR

The massage device cannot operate even though it has been connected to the power supply.

ADJUSTMENT

- Check the connection of the adapter.
- The device could be overheated after prolonged use. It will automatically cool down to resume its normal condition. Please turn off the power and allow the device to cool off for 30 to 50 minutes.

Care instructions

You have purchased a product manufactured with high quality materials. This requires appropriate care. For cleaning up the device just use a dry cloth, the silicone attachments can be rinsed under running water.

Storage

1. Recommended storage: dry, free of dust.
2. Keep liquids and open fire away from the Tappymed II.

Conformity declaration

Corresponds to the European guidelines EMC 2014/30/EU, LVD 2014/35/EU.

Warranty conditions

The dealer grants the statutory warranty period on the product described overleaf. The warranty applies from the date of purchase. The date of purchase is demonstrated by the sales receipt.

The dealer is to repair or replace products registered within Germany which are found to be defective free of charge. Parts subject to wear and tear such as cloth covers are not covered. To assert your claim, simply return the defective product with the sales receipt before the end of the warranty period.

The warranty claim does not apply if it is established that a defect has occurred e.g. through external influence or as a result of repair or alteration not carried out by the manufacturer or an authorised dealer.

The warranty granted by the seller is limited to repair or replacement of the product. Under this warranty, the manufacturer or seller has no further liability and is not responsible for damage caused by disregard of the operating instructions and/or improper use of the product.

If inspection of the product by the seller reveals that the warranty claim relates to a defect not covered by the warranty or that the warranty period has expired, the cost of the inspection and repair are to be borne by the customer.

Legal note: All contents do not contain healing statements. The diagnosis and treatment of illnesses and other physical disorders require a treatment by a physician, non-medical practitioner or therapist. The statements are strictly informative and may not be used as a substitute for a professional medical treatment. Every user is beholden to ascertain by accurate examination and, if required, after consultation of a specialist, if the application is the right one for their specific cases. Each dose or application takes place at the risk of the user. We dissociate ourselves from any healing statements or promises.

Casada International GmbH
Obermeiers Feld 3
33104 Paderborn, Germany
info@casada.com
www.casada.com

Copyright Images and texts are subject to the copyright of Casada International GmbH and may not be used for other purposes without express confirmation.

Copyright © 2017 Casada International GmbH. All rights reserved.

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Copyright © 2017 CASADA International GmbH.
Zdjęcia i teksty podlegają copyright firmy Casada International GmbH i nie mogą być wykorzystywane bez wyraźnej zgody.